



## QUICKTIP FLEXITUBE®



### THE CLEVER SOLUTION FOR CONSERVATION OF BOAR SEMEN



Millions and millions of sows have so far been inseminated with semen doses preserved in tubes made by Minitüb – the pioneer in artificial insemination technologies for sows. With Minitüb's many years of expertise in semen preservation, tubes provide a simple and effective solution within the semen processing line:

**semen production**



**transport**



**artificial insemination**

#### Your benefits

- easy handling
- easy to open by twisting tip without any tool
- suitable for all catheters and insemination techniques
- quick, easy and hygienic insemination
- smooth insemination: tube empties completely by itself – no semen waste
- optimum cooling rates and temperature stability at +17°C due to the well-balanced surface to volume ratio

**QuickTip Flexitube®**

REF.: 13452/xxxx






### Marking lines:

easily readable filling level

### Shape:

little air passages between the tubes allow for steady cooling of each semen dose in the fridge

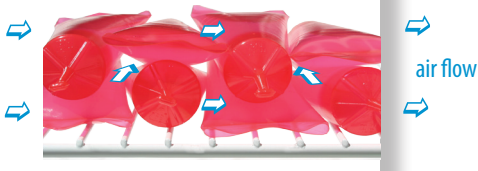

air flow

### Spout:

ensures an airtight fit of the tube on all different kinds of catheters




### Material:

sperm friendly polyethylene with superior cooling transfer ability

### Flexible walls:

for a complete semen uptake by the sow

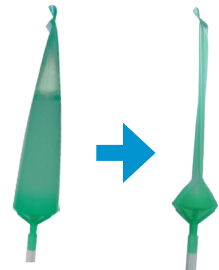

### Tube bottom:

particularly reinforced, so that the tube stands upright on the catheter

### Ergonomic twist-off tip:

no tool necessary with QuickTip tubes